

# SPIROEXCEL

PC based Pulmonary Function Test

Directly connect to USB port, no external power supply required



*digital turbine*

## Features

- Digital Spirometer with automatic Interpretation
- Easily sterilizable pneumotach
- Simple and Troublefree Technology
- Accurate and Precise
- Windows® based user friendly software
- Real time display of flow-volume loops
- Incentive screens for better efforts
- User's configurable 2 nos. prediction tables
- Pediatric Encouragement software
- Automatic BTPS Correction
- Multiple report formats with colour printout
- Unlimited test result storage and retrieval
- Report can be sent via e-mail



MEDICAID

An **ISO 9001 : 2000** Company



## Technical Specification

➤ Flow Measurement	: Bi-directional digital turbine	➤ Interface	: USB
➤ Volume Measurement	: Flow Integration	➤ Patients Incentive Screen	: Provided
➤ Flow Range	: 0.03-20 l/s	➤ Printer	: PC Compatible
➤ Volume Range	: 0-12 liters (BTSP)	➤ Storage	: On hard disk drive
➤ Flow Accuracy	: ± 2%	➤ Interpretation	: Automatic
➤ Volume Accuracy	: ± 2%	➤ Weight	: 100g.
➤ Graph Scale	: Meets ATs standards for graph	➤ Temperature	: 0-50°C
➤ Power Source	: 5V from computer through USB port.	➤ Dimension	: 165x50x25mm

### Prediction Table

Indian, ECCS, Knudson 83, Zapletal, Barcelona, ATS94, LAM, ITS White, ITS Black, LLN, user defined predicted calculations

### Test Performed

Forced Vital Capacity (FVC), Slow Vital Capacity (SVC), Maximum Voluntary Ventilation (MVV), Respiratory Pattern, Broncho-challenge tests.

## Measured Parameters




• FVC • FEV.5 • FEV 1 • FEV 3 • FEV1/FVC • FEV3/FVC • FEF 25-75 • PEF • FEF-25 • FEF-50 • FEF75 • SVC • ERV • IRV  
• TV • RR • PEFT • FMFT • FIV-1 • FIV1/FIVC • PIF • MVV

## Standard Accessories

➤ Spiroexcel Software CD	:01	➤ Carrying case	:01
➤ Operating Manual	:01	➤ Adapter for Pediatric mouthpieces	:01
➤ Digital Turbine	:01	➤ Adult mouthpieces	:10
➤ USB Cable	:01	➤ Pediatric mouthpieces	:10
➤ Nose Clips	:02		

## System Compatible Computer Hardware

Pentium-III, 128 MB RAM, 1.44 MB FDD, 40 GB HDD, CD-ROM, Keyboard, Mouse, USB Port, Colour Monitor 15",  
OPERATING SYSTEM : Window 98/ME/2000/XP. PRINTER : Any Inkjet Colour Printer supporting windows.

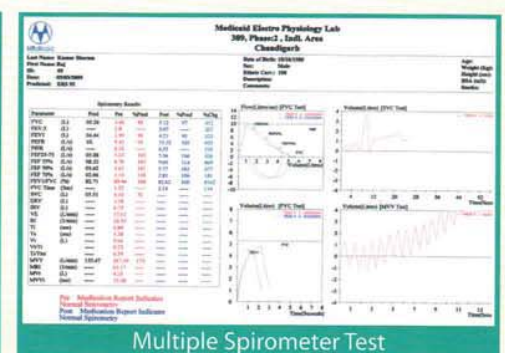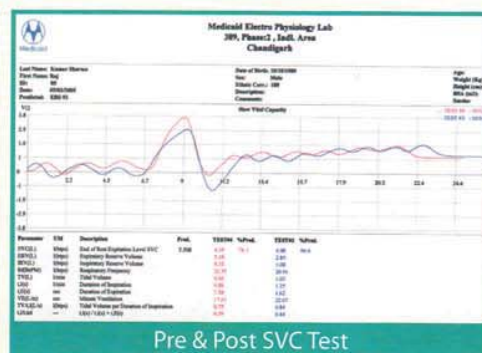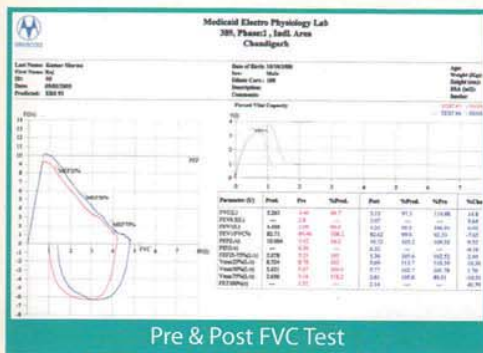

### Branch Office Ph No.:

Delhi	9810363242	Gwalior	9425112022
Ahmedabad	9824088842	Hyderabad	9849145151
Bangalore	9845191572	Chennai	9847230567
Lucknow	9839281804	Mumbai	9820929822
PB & HR	9814332314	Pune	9822619216
Jaipur	9414435131	Kerala	9840246486



## MEDICAID SYSTEMS

389, INDUSTRIAL AREA,  
PHASE-II, CHANDIGARH-160 002  
PHONE : 0172-2652706, 2641203 FAX : 2653608  
E-mail : medicaid@glide.net.in  
Website : www.medicaid.co.in